FINISHING TOUCH®
FLAWLESS
DERMAPLANE™

Features:


- 2 Long Precision Blades
- 2 Short Precision Blades
- Ridged Safety Guard
- Comfort Grip Handle
- 2 Caps

Flawless Dermaplane is to be used only ONCE a week.
Flawless Dermaplane should be used at night.

How to Use:

Flawless Dermaplane is recommended for exfoliation and hair removal on chin, forehead, cheeks and above the upper brow and upper lip area.

1. To begin, thoroughly wash and dry your face. The Flawless Dermaplane is for dry use only.
2. Select a clean Flawless Dermaplane for each use. Depending on the area, choose either the Long Precision Blade for larger areas, or the Short Precision Blade for smaller areas.
3. For best results, hold skin taut and blade at a 45° angle. Use short, feathery strokes and gently skim across skin. Do NOT press the Flawless Dermaplane down into the skin. No pressure should be used (See Images on Back). Only pass over the same area once during each use.
4. For tips and tricks and how-to-videos, please visit www.FlawlessProducts.ca

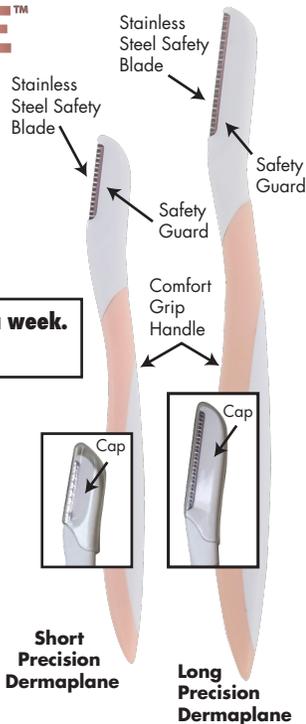After Use:

1. If skin is sensitive after use, cleanse your face with a gentle, non-acid cleanser.
2. Do not use lotion or make-up immediately after using your Flawless Dermaplane. Wait overnight before applying, to avoid clogging pores.

Warning:

- Do not use on skin that is experiencing any irritation including dryness, acne, sunburn, psoriasis, eczema, diabetes, circulatory or bleeding disorders.
- Avoid lips, nose, ears and eyes.
- Do not use if you have had a recent facial procedure.
- Not to be used for elimination of scars or tattoo removal.
- Never slide along the skin in the same direction as blade or it could cut the skin (See Images Above).

Distributed by
Mont-Royal, QC H4T 1X2
client@supertek.ca
www.clickit.ca
Customer service:
1-800-303-8354

FTDERMA4PK_090419ST

FINISHING TOUCH® FLAWLESS DERMAPLANE^{MC}

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