## INSTRUCTION MANUAL

- Gently removes stains, plaque, and tartar on contact.
- Cleans the front, back, and in between teeth.
- Promotes healthier gums.
- Recommended by Dental Professionals.

### To Install Battery:

- Locate the battery compartment cover on the bottom of the unit.
- Twist and pull the battery compartment cover off.
- Insert 1 "AAA" battery (not included) and replace the battery compartment cover by pushing and twisting until the two arrows are aligned.
- Do not mix old and new batteries. Do not mix alkaline, standard or rechargeable batteries


### To Use:

- Choose the attachment you'd like to use; either the stainless steel scaler tip, or the silicone flossing tip.
- To change attachments simply pull up on the head, remove, and insert the new attachment.
- Turn Sonic Pic on by pressing the button on the side, and then clean and whiten your teeth by simply pointing and pressing!

### To Clean:

#### **IMPORTANT: Do not submerge in water.**

- Make sure Sonic Pic is turned off before attempting to clean.
- Clean the attachments by removing them from the unit and cleaning them in water or alcohol pads.

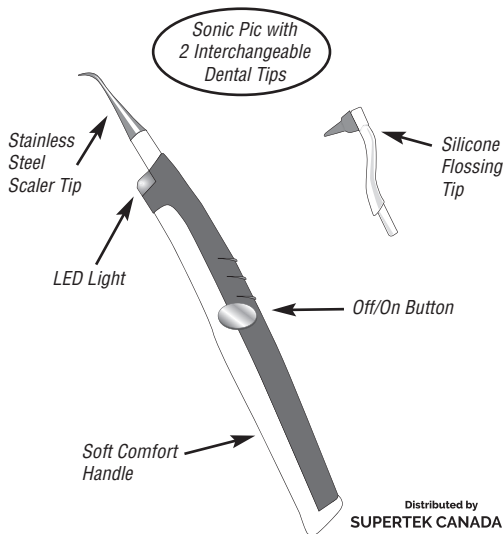Made in China

Distributed by  
SUPERTEK CANADA Inc.  
Mont-Royal, QC H4T 1X2

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